

Artificial eye technique





Mizushima prosthetic eye lab.
Fumio Mizushima (JAPAN)
mizu-eye@wing.ocn.ne.jp

Pigmented iris button making

Making artificial eyes

Hand made buff and polishing

Making iris button

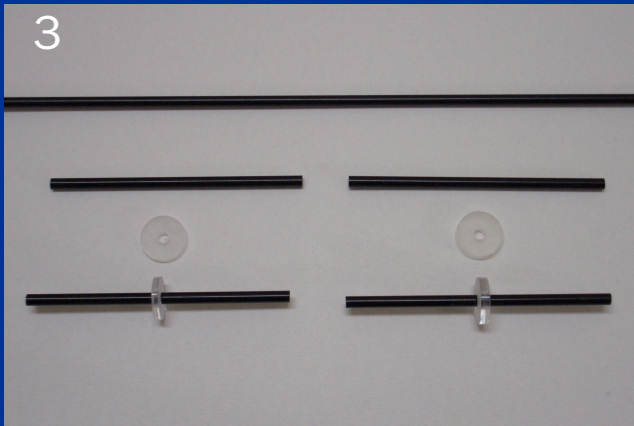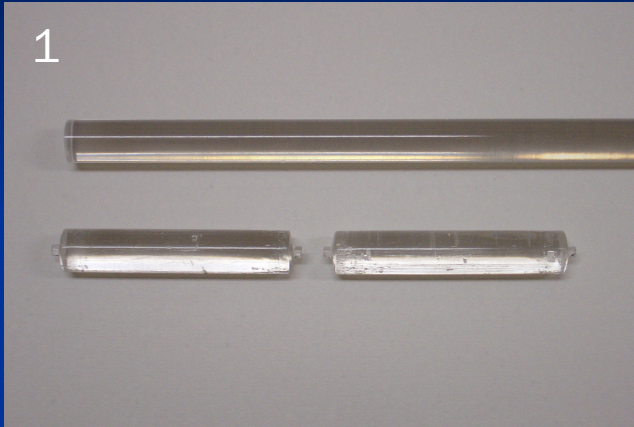

Making iris button



mixing materials

Making iris button



First curing

Making iris button



Processing with a lathe

Making iris button



Second curing

Making iris button



Slice with a lathe

Stock variety iris buttons




Pigmented iris button making

Making artificial eyes

Hand made buff and polishing

Processing iris button to get natural limbus

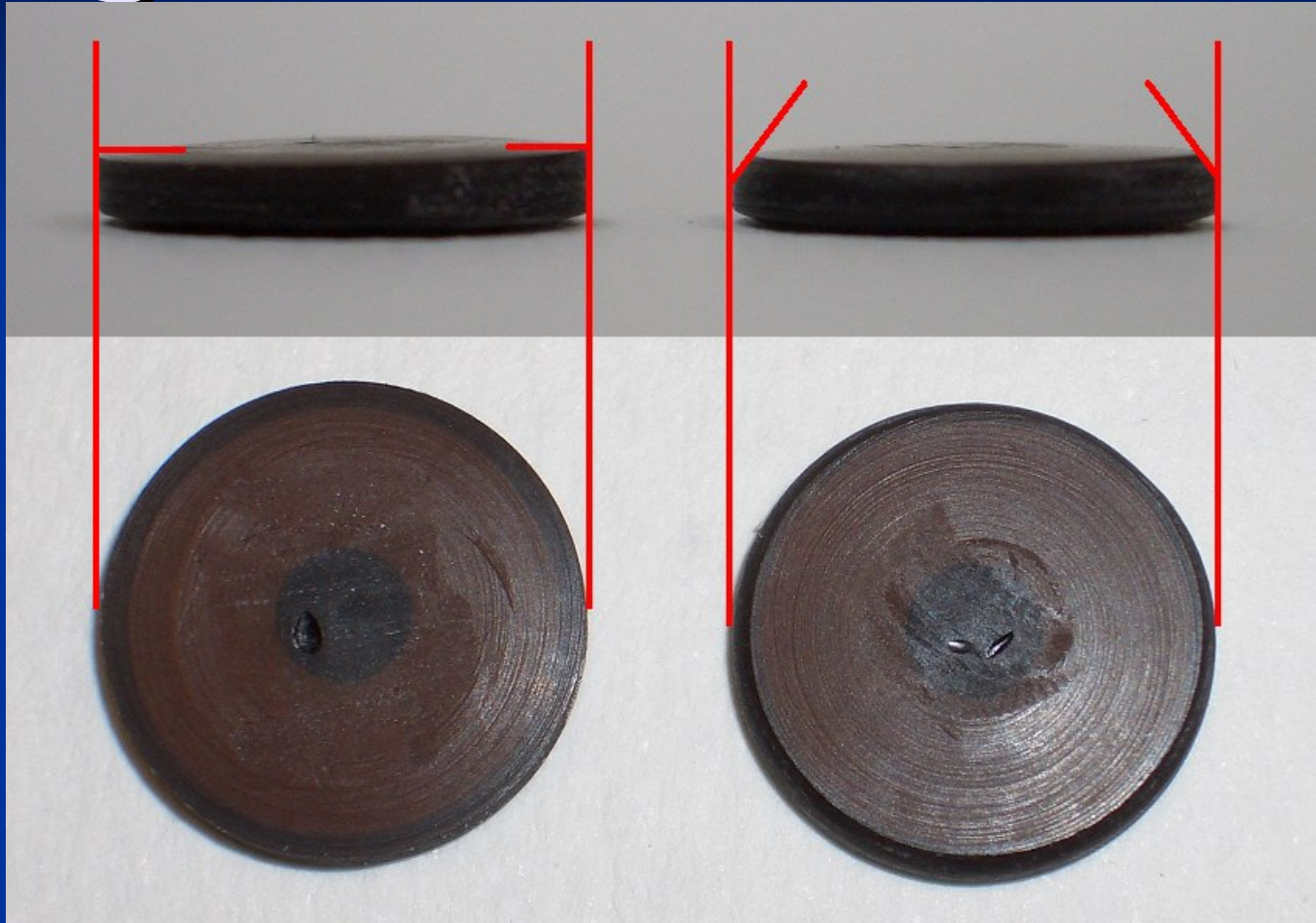

Remove iris button edge



press iris button to curve in boiled water

Fixing iris button to form



Press the form





Process of heat to curing




A pressure cooker was used for curing

Grinding sclera layer

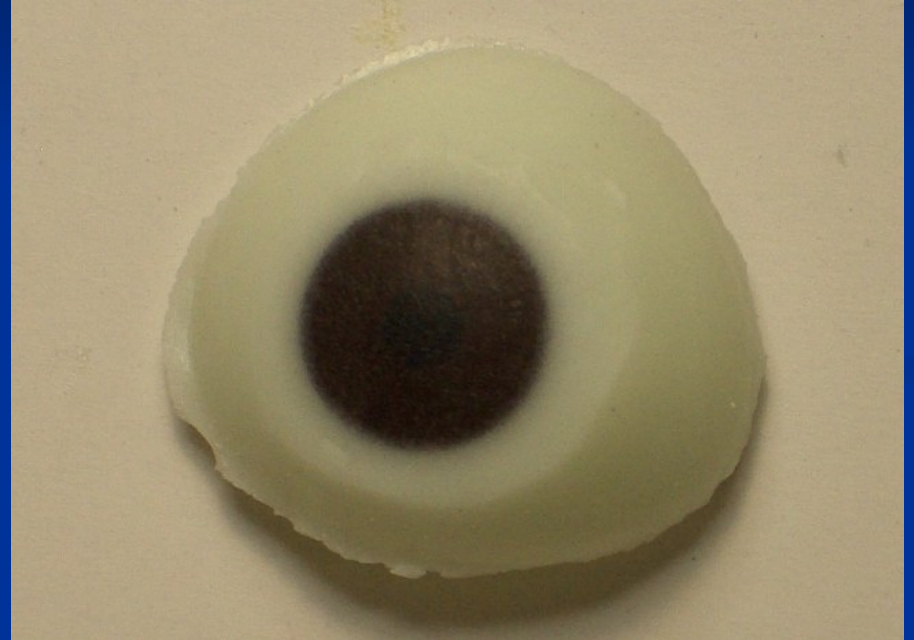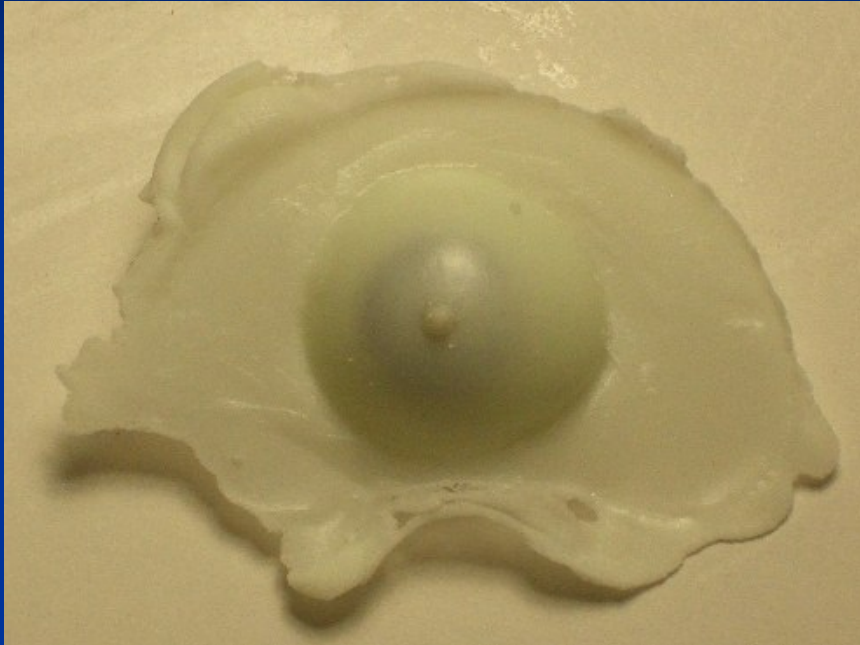

Process making iris natural


1. Make lines on iris
by scratching

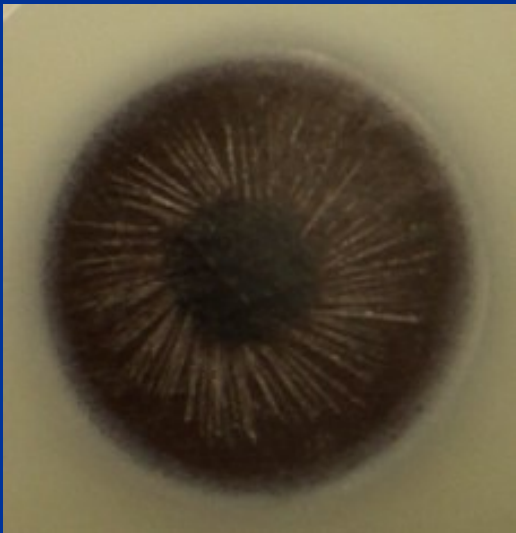


3. Make the surface bumpy
to achieve reality



2. Fill black ink in lines

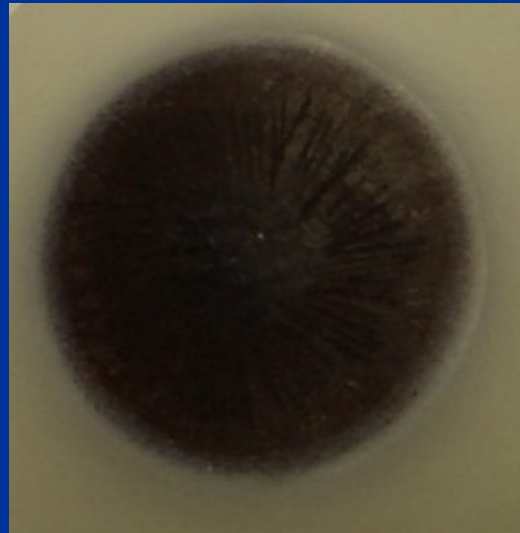

Making blood vessels and covering with clear acrylic materials

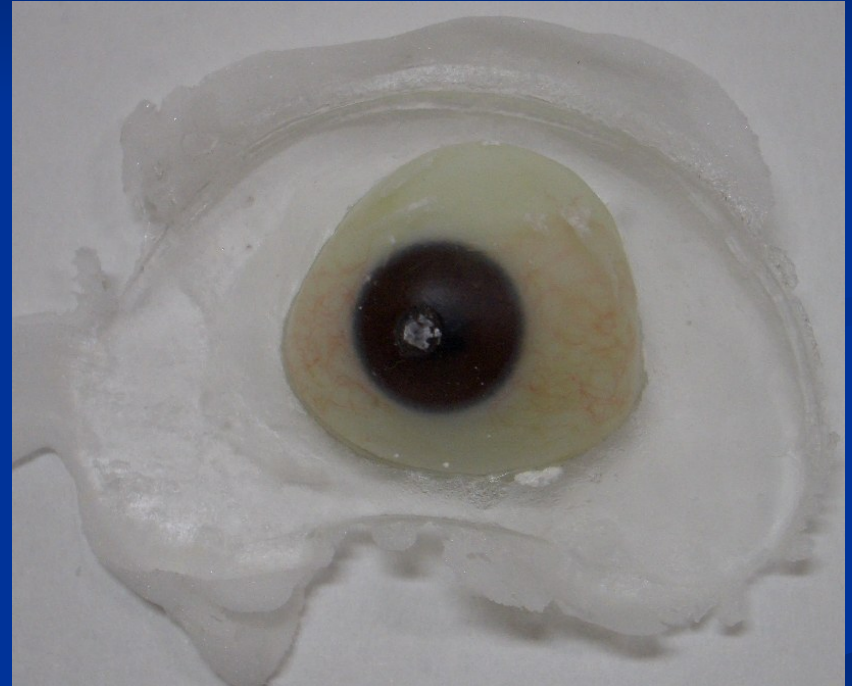

Pigmented iris button making



Making artificial eyes

Hand made buff and polishing

Polishing tools




Hand made rough polisher

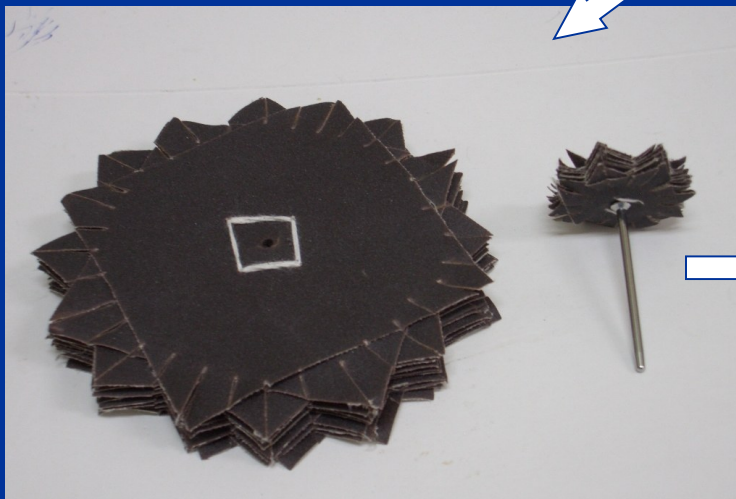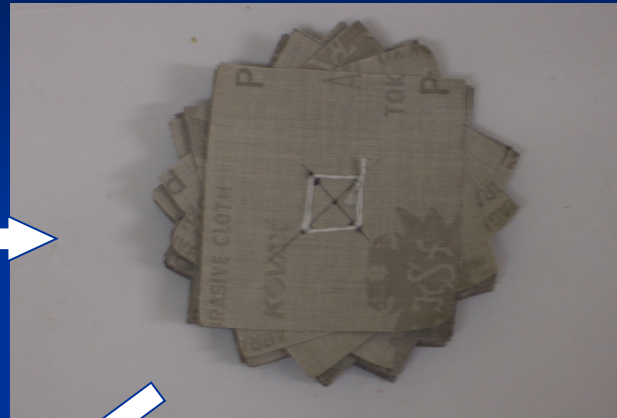

Middle polisher




Final polisher

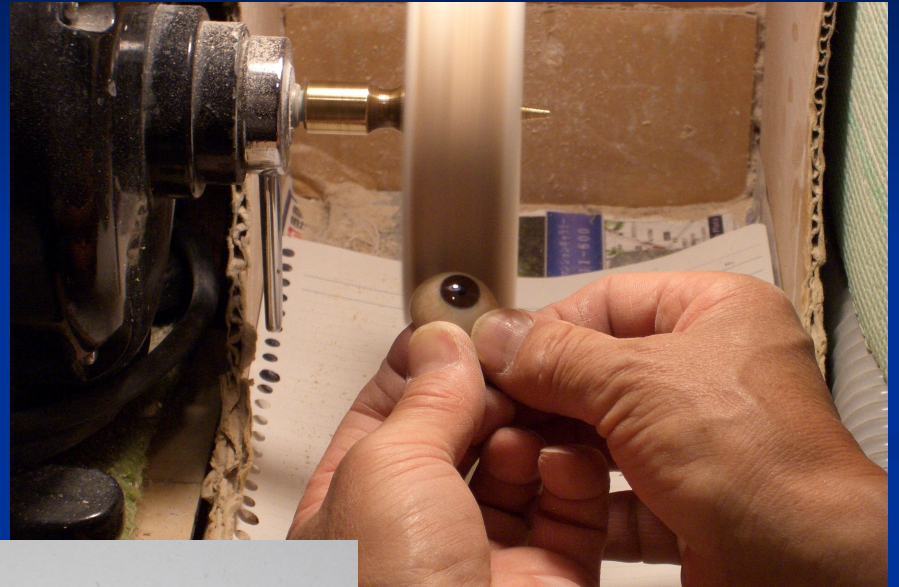




Wearing artificial eye

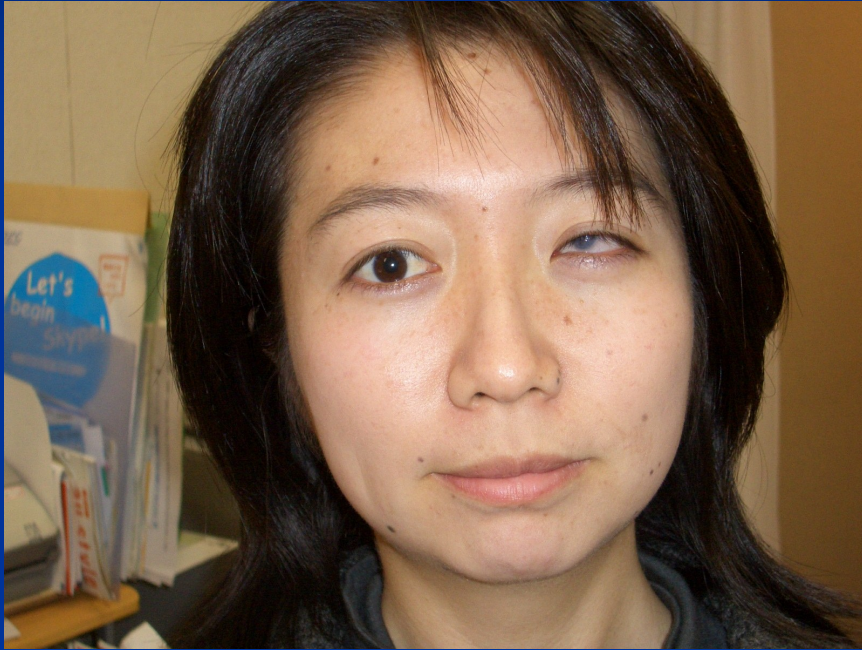


Before

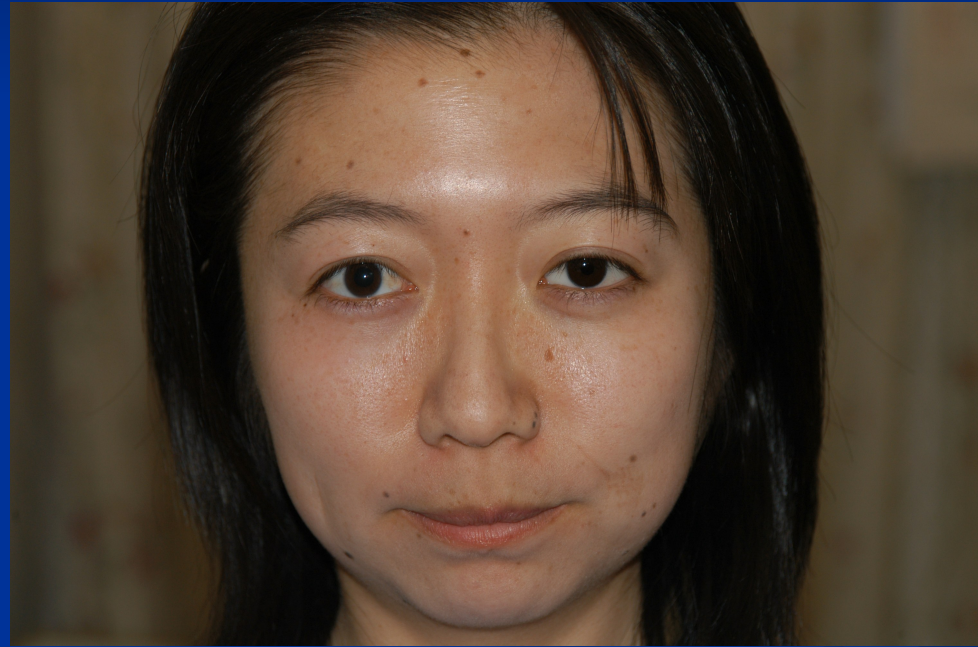

After

Wearing artificial eye

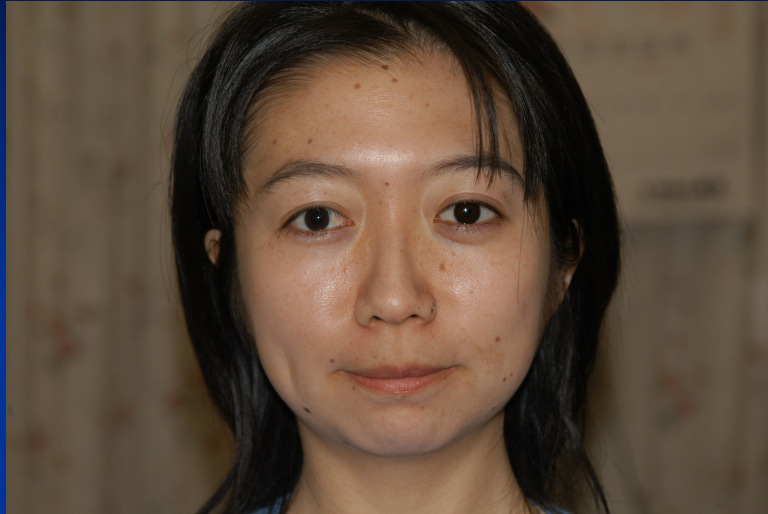

Wearing artificial eye

